


## Instruction Manual

Emjoi®  
**MICRO-Pedi™**  
*Power*



AP-3RPA  
Electrically-operated

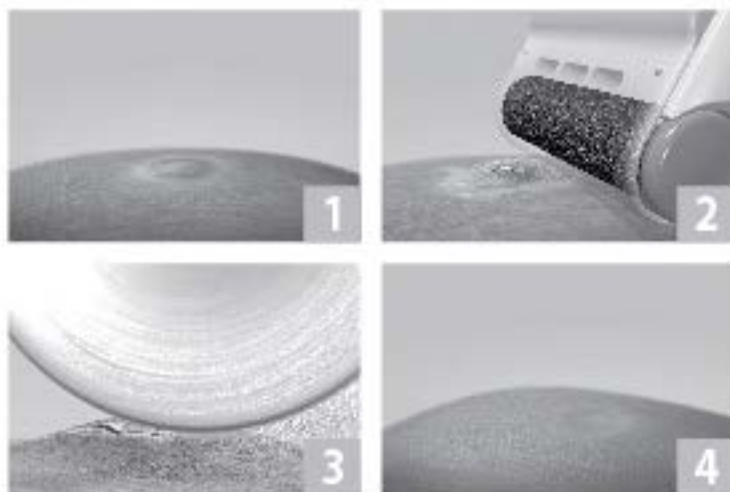

# GENERAL DESCRIPTION OF PARTS



1. Protective Cover
2. Micro-Pedi Power Roller
  - blue for regular skin
3. End Cap of Micro-Pedi Power Roller
4. Safety Switch Lock
5. On/Off switch
6. Release Button for Micro-Pedi Power Roller
7. Adapter Inlet



The **Emjoi® Micro-Pedi Power** is a super powerful, electrically-operated device created to smooth dry, rough and calloused skin on the feet. Unique roller spins 360°, at an amazing 30 times a second, to gently and effectively buff away hard, jagged skin without the hassle of constant rubbing. The **Micro-Pedi Power** gives salon results in seconds.

**Warning:**

The unit is very powerful. Gently rotate (back and forth or side to side) over the callous / corn / hard skin area for 2-3 seconds. Stop and check if you achieved the softness you desire. If not, go over the area for 2-3 more seconds and check again. Continue to use this method until the level of smoothness is obtained. Over use may result in bleeding or injury.

**Caution:**

Do not store or use your **Emjoi® Micro-Pedi Power** where it can fall into a tub or sink of water. Do not reach for your unit if it has fallen into water, immediately unplug it.

Each **Emjoi® Micro-Pedi Power** roller is recommended for one person only. For hygienic reasons, the **Micro-Pedi Power** roller should not be shared with others.

# INSTRUCTIONS FOR USE

1. Insert the adapter plug into the adapter inlet and plug your adapter into the electrical outlet (fig. 1 & 2).
2. Remove the plastic protective cap (fig. 3).

**Caution:**

Before plugging into an electrical outlet, make sure that the adapter is appropriate for the local voltage. This is extremely important when traveling overseas.



# HOW TO USE

Emjoi® Micro-Pedi Power is designed for smoothing coarse skin on the feet and heels (fig. 4). Do not use this device anywhere else on the body.

**Note:** Only use Micro-Pedi Power roller on dry skin.

1. Wash and clean the skin and dry fully.

**Note:** Do not use on broken or bleeding skin.

2. Make sure the Micro-Pedi Power roller is placed into the unit. See section of **Changing the Micro-Pedi Power Roller** on how to remove and replace rollers.
3. Remove the protective cover from the unit and it is ready to use.



4. Turn on the unit by pressing the safety switch lock, in the center of the and at the same time push the switch up (fig. 5).

**Warning:** The unit is very powerful. Gently rotate (back and forth or side to side) over the callous/corn/hard skin area for 2-3 seconds. Do not press hard - gently glide the unit over desired area. The unit will stop if it is pressed too hard on the skin. Stop and check if you have achieved the softness you desire. If not, go over the area for 2-3 more seconds and check again. Continue to use this method until the level of smoothness is obtained. Over use may result in bleeding or injury.

**Warning:** Stop using if skin becomes sore or inflamed.

**Warning:** Do not use on skin for more than 2-3 seconds at a time.

5. Turn the unit off (fig 6).
6. Rinse the skin or use a wet towel to get rid of the excess dead particles.
7. To clean the **Emjoi® Micro-Pedi Power**, please see the Cleaning section.
8. Keep the protective cover on your unit when it is not in use.





# CHANGING THE MICRO-PEDI POWER ROLLER

Change the Micro-Pedi Power roller as needed for best result. Roller last for a few sessions depending on frequency of use and size of area being treated. Replace as needed.

**Replacement Emjoi® Micro-Pedi Power rollers are available at most retailers that sell the Emjoi® Micro-Pedi Power unit.**

The Micro-Pedi Power roller contains fragile parts, please handle them with care.



Please take extra care when replacing the Micro-Pedi Power roller.

**Note:** Micro-Pedi Power rollers that are broken or worn could cause damage to the skin.

Replace the **Emjoi® Micro-Pedi Power** roller when you find any defects to the roller or when you find them less effective.


Make sure the unit is turned off before replacing the Micro-Pedi Power roller.

1. Pull down the release button, on the side of the unit, while the end cap is facing up and the Micro-Pedi Power roller will pop out (fig. 7 & 8).
2. After the Micro-Pedi Power roller pops out, remove it from the cylinder holding the end cap of the Micro-Pedi Power roller (fig. 9) and dispose of it.

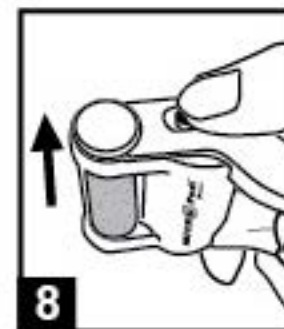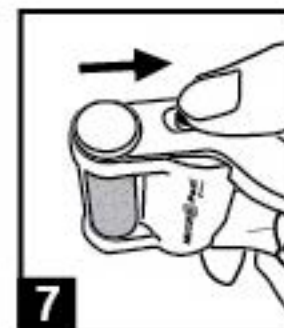

3. Remove the Micro-Pedi Power roller from the replacement container (fig. 10), (replacement rollers sold separately), by pressing the bottom with your pointer finger and catching it in the palm of your opposite hand.
4. Place a new Micro-Pedi Power roller into the unit until you hear a 'click' sound. Be sure to insert the new Micro-Pedi Power roller into the unit with the thinner side going into the unit (fig. 11 & 12).
5. After placing a new Micro-Pedi Power roller into the unit check that it has been firmly installed with no signs of shaky or loose parts (fig. 13).

**Micro-Pedi rollers can also be used with the Micro-Pedi Power.**

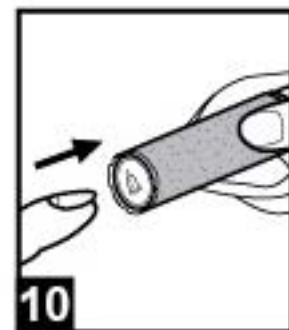



# CLEANING

The **Emjoi® Micro-Pedi Power** should be cleaned after each use. It is not water resistant so do not submerge or rinse the unit in water as it will cause defects to the unit. **Micro-Pedi Power rollers should not be placed in water.**

1. Before cleaning, make sure to turn off your unit and unplug it from the electrical outlet.
2. Remove the protective cover from the unit.
3. Use the brush to remove any debris from the head of the unit.
4. Wipe the body with a clean dry cloth.
5. Replace the protective cover.

# WARNING / IMPORTANT

- The unit is very powerful. Gently rotate (back and forth or side to side) over the callous / corn / hard skin area for 2-3 seconds. Stop and check if you achieved the softness you desire. If not, go over the area for 2-3 more seconds and check again. Continue to use this method until the level of smoothness is obtained. Over use may result in bleeding or injury.
- Before using your **Emjoi® Micro-Pedi Power** please read all these instructions carefully and completely.
- The **Emjoi® Micro-Pedi Power** is ONLY for use on the feet and heels. **Do not use this device anywhere else on the body.**
- Use the **Emjoi® Micro-Pedi Power** only for purposes described in this instruction book.
- Please keep these instructions for future reference.
- This unit is intended for one person. Do not share this unit with others.
- Each **Emjoi® Micro-Pedi Power** roller is recommended for one person only. For hygienic reason, the Micro-Pedi Power roller should not be shared with others.
- Do not use if you have diabetes or poor blood circulation.
- If you have any skin conditions or other medical issues, please consult with your physician prior to use.
- Do not use on skin for more than **2-3** seconds at a time.
- Do not use the Micro-Pedi Power roller when you find worn or nicked edges on the roller.
- Change the Micro-Pedi Power roller as needed for best result. Roller last for a few sessions depending on frequency of use and size of area being treated.
- Do not use on skin that is red, irritated, inflamed, infected or blistered.
- Do not use on warts.
- Do not use on sores or open wounds.
- Discontinue use immediately if you experience any pain or irritation.
- Overuse may cause skin irritation.
- Seek medical attention if skin becomes severely irritated or bleeds.
- Only use the **Emjoi® Micro-Pedi Power** on dry skin.
- This unit is for dry use only.
- Household use only.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety.
- Do not place or store your unit in direct sunlight. Please store your unit in a cool dry place at a temperature between 0 - 40°C.
- Keep appliance dry. It is **not water resistant** so do not submerge or rinse the unit in water as it will cause defects to the unit.
- Do not store your unit where it can fall into a tub or sink of water.
- Never use your unit if it has been damaged or submerged into water.
- Do not use water or any other liquid on the Micro-Pedi Power unit or roller.



- Do not turn your unit On if the Micro-Pedi Power roller is not attached.
- Always turn your unit "Off" when it is not in use and before cleaning it.
- Do not use this unit if the Micro-Pedi Power roller is damaged as this could lead to injury.
- Do not drop or insert anything into any openings of the unit.
- Do not use your unit outdoors.
- Do not use your unit where aerosol products or oxygen is being used.
- Children should be supervised to ensure that they do not play with this appliance.
- Keep your unit out of reach of children as it is not intended for use by children.
- This appliance is not to be used on, by or near children or invalids.
- To prevent accidents and injury to you or damage to the unit – when your **Emjoi® Micro-Pedi Power** is in use, keep it away from the hair on your scalp, your eyebrows, eyelashes as well as clothes, brushes, wires, cords, shoelaces, etc.
- If you have any medical conditions, issues, questions or any sensitivity whatsoever we strongly advise you to consult with your health care provider before using this product.
- If you have a pacemaker consult with your health care provider prior to use.
- Do not continually rub the unit on the same area or apply strong pressure.
- Keep the protective cover on your unit when it is not in use.
- Place the adapter into the unit and then into the electrical outlet. Make sure that the plug is inserted firmly into the appliance. When you are finished using your unit turn the unit "OFF" and then remove the adapter from the electrical outlet and then from the unit.
- Do not reach for your unit if it has fallen into water. Immediately unplug the unit.
- When you unplug your unit please pull on the plug and not on the cord.
- Always unplug your unit when it is not in use and before cleaning it.
- Do not wrap the cord of your adapter around the unit.
- Keep the cord of the adapter away from hot surfaces.
- Never use your unit if the adapter cord or electrical outlet is damaged.
- Do not cut off the adapter to use it with another plug as this can be a hazardous situation.
- Only use your Emjoi® Micro-Pedi Power with the adapter supplied with the unit. For questions regarding your adapter in the US contact 1-888-99EMJOI, M-F, 9-5 EST.
- Caution: Before plugging into an electrical outlet, make sure that the adapter is appropriate for the local voltage. This is extremely important when traveling overseas.

## Warranty

The **Emjoi® Micro-Pedi Power** warranty against defects in materials and workmanship lasts one year from the date of purchase.

**Keep original receipt for warranty service.**

During this warranty period, return the product to Emjoi, Inc. with proof of purchase to FDR Station, PO box 7406, New York, NY 10150. Emjoi will take care of repair, or at their discretion, replacement of a defective unit without charge. For questions regarding this policy you may call toll free 1-888-99-EMJOI (1-888-993-6564) M-F, 9-5 EST.

**Important:** This warranty does not cover damages resulting from accidents, misuse, alterations, abuse, lack of reasonable care, or performance of services or repairs by an unauthorized service center.

The manufacturers and distributors of the **Emjoi® Micro-Pedi Power** do not assume responsibility for any damage caused by misuse or use other than that described in these instructions.

Model: AP-3RPA

Patents pending

Made in P.R.C.

© Soft Lines International Ltd. 2013

**w w w . e m j o i . c o m**

PK-IB-EMJ-3RPA-01-0